



Inbios International, Inc.
562 1st Ave. South
Suite 600
Seattle, WA 98104
206-344-5821 (phone)
206-344-5823 (fax)
www.inbios.com

Product Specification Sheet

Product	Monoclonal antibody specific to Zika Virus NS1
Catalog Number	INB9162
Description	Antibody specific to Zika Virus NS1
Antigen	Recombinant Zika Virus NS1
Clone	6613.7A1
Host Animal	Mouse
Format	Purified, liquid
Isotype	IgG1
Buffer	PBS pH 7.4
Application	Suitable for use in ELISA (Fig 1). Each laboratory should determine an optimum working titer for their own application. Other formats have not been tested.
Storage and Handling	Store at -20°C. Avoid multiple freeze/thaw cycles. Prepare working dilution only prior to immediate use. Centrifuge before opening to ensure complete recovery of vial contents.

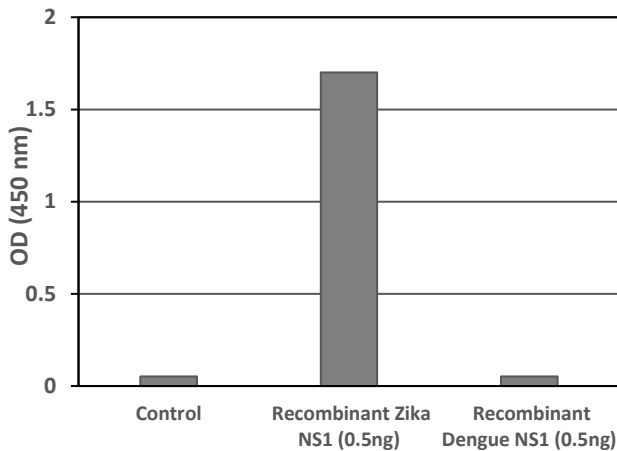


Fig 1. Preliminary ELISA performance 1µg/mL Mouse monoclonal antibody to Zika Virus NS1 (Clone 6613.7A1) on normal human serum (control) and normal human serum spiked with either 0.5ng recombinant Zika virus NS1 or Dengue virus NS1.

For Research Use Only. Not for use in Diagnostic Procedures.